



2,000~3,000 shot

2 times

4~6 treatment

No downtime

One session

A week

Total treatment session

TECHNICAL SPECIFICATION

Air pressure(bar)	1~5 bar (0.1step)
Frequency	1~16Hz (1step)
Head	15mm, 20mm, 35mm
Dimension(mm)	W385 x D450 x H980
Weight	40kg



Rm.205~207, Woolim e-biz center II, 12, Digital-ro 33-gil, Guro-gu, Seoul, Korea [08377]
 TEL. 02-871-0657 FAX. 02-871-0658 URL. www.k1med.co.kr E-MAIL. sand812@k1med.com

EXTRACORPOREAL SHOCKWAVE
THERAPY


ComboWAVE

ESWT
DUAL RADIAL TYPE



EXTRACORPOREAL SHOCKWAVE THERAPY

Weighty bullet is accelerated by air pressure to transmit kinetic energy to the head.
And kinetic energy is converted to shockwave in the head and spread radially from the contact area with the skin inside the body.

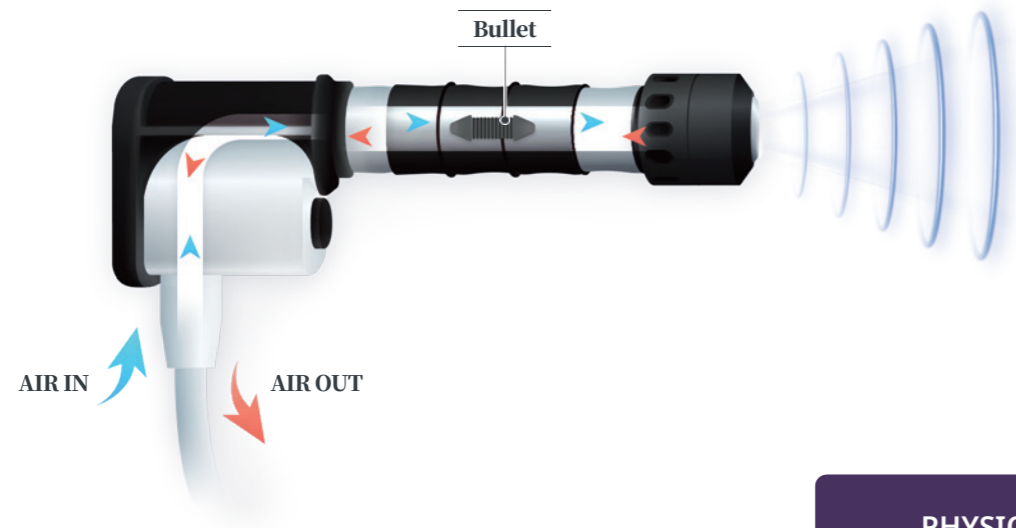

PHYSIOTHERAPY

Calcification breaking
 Tendons, Muscle pain relieving
 Revascularization promoting
 Calcareous removing
 Prostatitis care



PHYSIOTHERAPY

COMBOWAVE treatment is a non-invasive treatment. Involves transmitting high-energy sound waves or sound energy to the affected area of the body, thereby triggering the body's own natural repair mechanisms and simulating healing.

ESWT WORKS INSIDE BODY

- Stimulate fibroblast and tenocyte production
- Prompt Body Response : increase blood circulation
- Increase metabolism in the injured tissue
- Activate cell generation
- Dissolve calcium deposits
- Eliminate causes of chronic pain in a short period of time



COMBOWAVE INDICATION

